


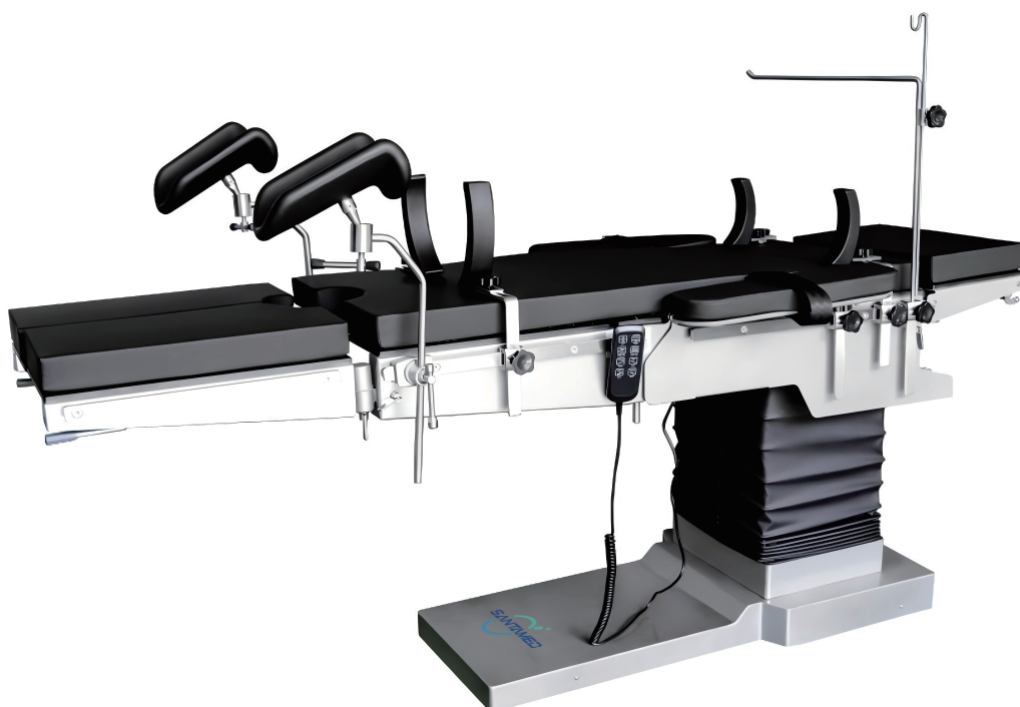

STable-8600

OT/OR Table

- Unique single-side column design with ultra-thin base, which reserves extra space for C-arm, U-arm, X-ray fluoroscopy and other equipment.
- Decompression mattress with antistatic, waterproof, and seamless design removable.
- 8-direction position movement with anti-clash function by remote control.
- Head and leg board can be exchanged according to different surgery demands.
- Return to original position with one button optional .
- Foot treadle break device easy for moving and operation.


Features:

- Mattress is wrapped with quality leather, high-density memory sponge in one piece, with seamless and antistatic, waterproof design.
- Various actions are controlled by low voltage DC motor with advanced backup power, which is available to move the table body max 300 mm in totally 8 directions.
- MCU processing is used, self-locking function is provided to avoid false triggering.
- Quality stainless steel is used for table support, base plate, column shield and rails on both sides of the table, easy for cleaning and sterilization.
- Table support system is with heavy weight capacity up to 200 kg.
- Famous and quality carbon steel plate is used for table top, which can be used for C-arm/X-ray checkups and radiology imaging.
- Special design of foot treadle brake device easy for moving and fixed position operation use.
- Abundant accessories applicable for multiple position use applicable for different operations.

ST Series professional operation table parameters and specifications:

STable-8600			
Table Length (mm)	2135 ± 20	Backplane Downward	≥ 15 °
Table Width (mm)	550 ± 20	Backplane Upward	90°
Height Adjustment (mm)	625-925 ± 20	Head Board Fold Upward	≥45°
Loading Capacity (kg)	200	Head Board Fold Downward	≥80°, detachable
Leg Turning Angle Upward	≥15°	Lumbar Rise Height (mm)	≥ 130±10
Leg Turning Angle Downward	≥ 60°, removable		
Leg Turning Angle Fold Out	≥ 90°	Angle of the Table Top	≥ 20 °
Backward Tilt Forward	≥ 21 °	Transparent Range	≥ 1350
Backward Tilt Reclining	≥ 21 °	Main Voltage (AC)	220V 50Hz

Accessories	Quantity
Anesthesia Screen	1 item
Shoulder Support	2 pcs
Foot Board	2 pcs
Screen Frame Slider	2 pcs
Armrest	2 pcs
Slider	6 pcs
Support Frame	2 pcs
Leg Support	2 pcs

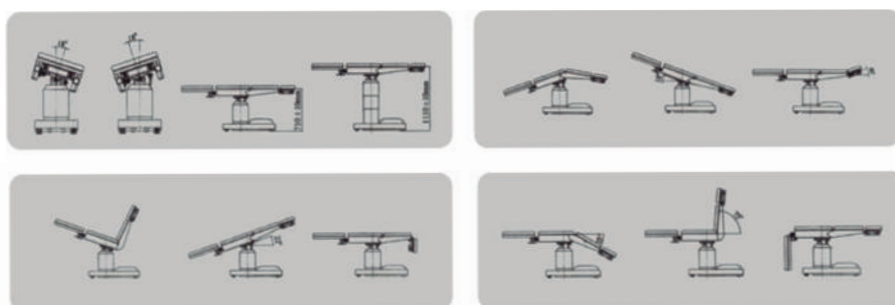

SANTAMED

ZHENGZHOU SANTA MEDICAL CO.,LTD

🏠 Room 301 Building-6, Jincheng Time Square, Jinshui Zone, Zhengzhou450000 Henan Province PR,China

☎ + 86-177 19820906

✉ sales@santamedhk.com; cloverye@santamedhk.com

www.santamedhk.com